


Gerber Instruments

NEW – Testing whipped cream with the Gerber CT III

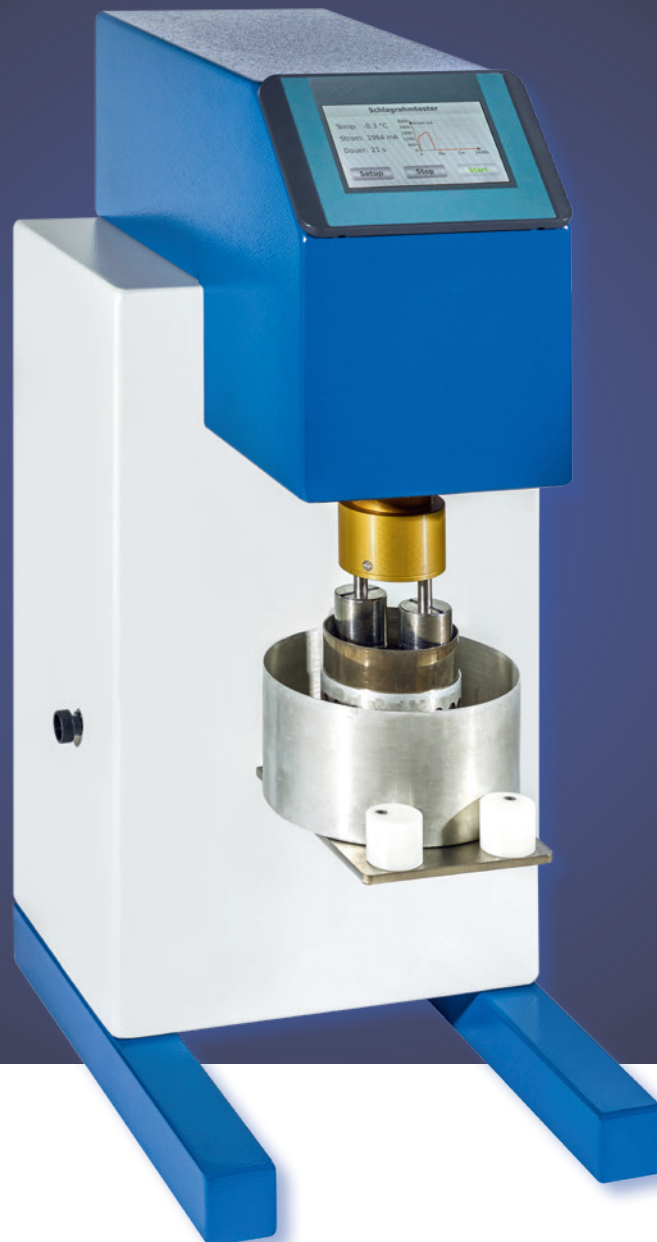

GERBER
INSTRUMENTS
Since 1892 

Testing whipped cream with the Gerber CT III

Since its foundation in 1892, the company Gerber Instruments manufactures devices and analytical technology for the quality control of the dairy industry.

The testing device for whipped cream CT III, which has been newly developed by the company Gerber Instruments, supports the daily routine controls in the laboratory and production development. It delivers valuable measurement results for the safety and quality of your products. With the new Gerber CT III the physical test of the cream that needs to be whipped will be performed with the following properties:

- Whipping capability
- Maximum attainable whipping resistance
- Optimum whipping time
- Formability and settling (*with optional accessories)



- Measuring method of the CT III made by Gerber
The testing instrument for whipped cream CT III works according to Prof. Dr. Mohr, the DLG as well as the Swiss Food Code. The cream sample will be whipped exactly under defined conditions: standardised sample beaker, standardised whip, constant speed
- Novelties of the CT III made by Gerber
The measuring instrument is operated by a touch screen. With the help of the graphic on the display you can track the course of the measurement on the display. The drive is equipped with a digital speed control. The sensitivity can be set as required. It is possible to output the data on the LIMS or a printer.

Technical data

Sample quantity	100 ml
Sample temperature	5 °C
Display of the torque	in mA or Nm
Display of the whipping time	in seconds
Mains connection	115 / 230 VAC, 50 / 60 Hz
Dimensions (W × D × H)	270 × 355 × 630 mm
Weight	19 kg
Guarantee	2 years

Order No.	Type
09.1002	Tester for whipped cream CT III incl. sample beaker with cam (nickel-plated), cooling jacket (aluminium) as well as 2 whips (stainless steel)
09.2010*	Settling mould plate to mould a whipped cream tube of 6 × 6 cm for the settling sample
09.2011*	Settling screen 10 × 10 cm according to DIN to determine the settled amount of liquid
09.2020*	Dipstick to determine the increase in volume of the whipped cream, scale: from 100 to 250 ml

Gerber Instruments AG

Im Langhag 12 - CH-8307 Effretikon (ZH)
info@gerber-instruments.ch - www.gerber-instruments.com
Phone +41 (0)52 343 37 37 - Fax +41 (0)52 343 30 70

GERBER
INSTRUMENTS
Since 1892 